

PCR protocol

RBRC 02219

Strain name B6;129-Use1 <tm1.1>

Primer	primer_name	primer_seq	length	
1	p31 geo F	TTAAACCTCAGCACTGAGGAGGCAG	25	Depositor's primer
2	p31 WT R	AGGCAAAGATGACCCTGGGACCTCT	25	Depositor's primer


Reaction Components

Enzyme	TaKaRa Ex Taq (Cat. No. RR001A)
--------	---------------------------------

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		17.3			
PCR buffer	10 X	2.5			
dNTP	2.5 mM	2.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
DNA		1.0			
Taq		0.2			
total		25			

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
3	63	30*	
4	72	60*	*Repeat for 30 cycles
5	72	300	

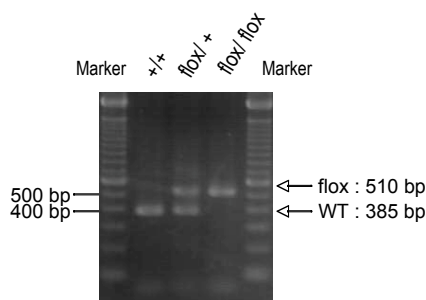

Marker

100 bp ladder (Invitrogen, Cat. 15628-050)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set		product_size	product
primer 1	primer 2	385 bp	WT
primer 1	primer 2	510 bp	flox

Comment