

PCR protocol

Sheet ID PS_03262

Release date : 08/01/2012

Print date : 8/1/2012

RBRC 04180

Strain name C57BL/6N-Tg(Thy1-ChR2/GFP)437-5Sasa/Rbrc

Primer	primer_name	primer_seq	length	
1		GCTGGAATCAAAGGTGTGAGCTC	23	RIKEN BRC
2		CTTGTACAGCTCGTCCATGCCGAG	24	RIKEN BRC
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol


Reaction Components

Enzyme TaKaRa LA Taq with GC buffer (Cat. No. RR02AG)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		3.3			
PCR buffer	2 X GC I	10.0			
dNTP	2.5 mM	1.6			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
primer 4	10 pmol/ μ l	1.0			
DNA		1.0			
Taq		0.1			
total		20			

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
3	60	30*	
4	72	120*	*Repeat for 30 cycles
5	72	300	

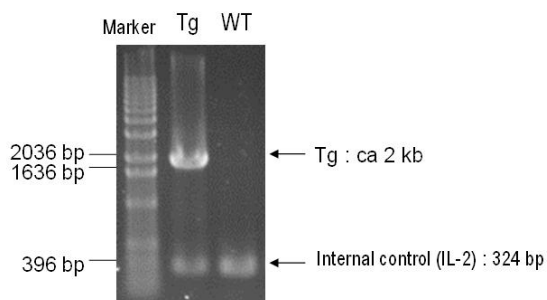

Marker

1 kb DNA ladder (Invitrogen, Cat. 15615-024)

Gel

1% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	primer 2	ca 2 kb
primer 3	primer 4	324 bp

Comment