

PCR protocol

Sheet ID PS_06801
Release date : 7/14/2016
Print date : 7/14/2016

RBRC 04909

Strain name B6;129P2-Il1f6<tm1Yiw>

Primer	primer_name	primer_seq	length	
1	IL1F6-common	AGGTTAGAAACGCAGCCTGT	20	Depositor's primer
2	IL1F6-wild	TTGCTCTTTCCTTGGGACTG	20	Depositor's primer
3	IL1F6-mutant	CGGCCATGATATAGACGTTG	20	Depositor's primer

Reaction Components


Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		1.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Standard PCR

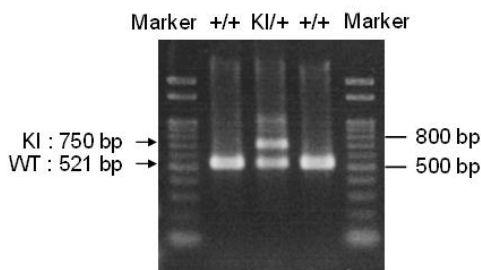

Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set		product_size	product
primer 1	primer 2	521 bp	WT
primer 1	primer 3	750 bp	KI

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.