

RBRC 05661

Strain name B6.Cg-Pcp4<tm1.1>

[illegible]


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

[illegible][illegible]

Standard PCR

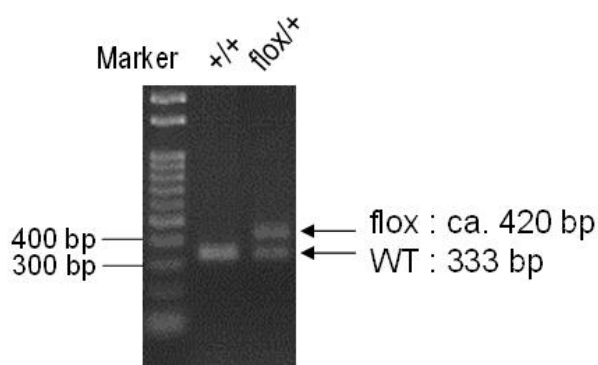

Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set		product_size	product
primer 1	primer 2	333 bp	WT
primer 1	primer 2	ca. 420 bp	flox
primer 2	primer 3	(570 bp)	null

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.