

PCR protocol

Sheet ID PS_07239
Release date : 5/17/2017
Print date : 5/11/2017

RBRC 09923

Strain name C57BL/6N-Zfhx1a<tm3.1Yhi>

Primer	primer_name	primer_seq	length	
1	intron5 MN fwd	GCAGATCCATGCATATTGATAGGCTGATGG	30	Depositor's primer
2	intron5 B6 rev	GTGATGAATATGGCTCTGCACCTCGGACAG	30	Depositor's primer
3	intron6 B9 rev	AGAACCTCAAGATGGTCTCAGAGAAGTCAC	30	Depositor's primer

Reaction Components


Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
--------	-----------------------------------

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		1.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
primer 3	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.2			
total		20			

Standard PCR

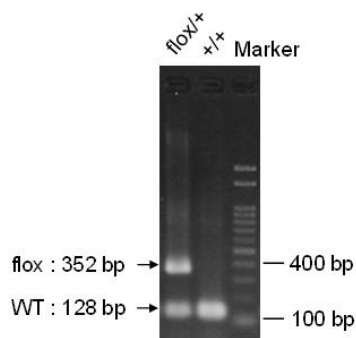

Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1	128 bp	WT
primer 1	352 bp	flox
primer 1	276 bp	null (After cre mediated excision)

Comment

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.