



Sheet ID PS_07666
Release date : 11/6/2018
Print date : 11/6/2018


Strain name C57BL/6-Slc6a5^{<tm1.1(cre)Ksak>}

[illegible]

Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
---------------	-----------------------------------

[illegible][illegible]

Marker	Gene Ladder 100 (Wako, Cat. 316-06951)
Gel	2% agarose in 1 X TAE (100 V, 30 min)

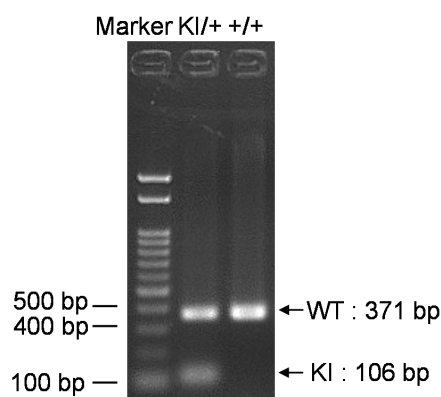

primer_set		product_size	product (allele type)	
primer 1	primer 2	371 bp	WT	+/+ : 371 bp KI/+ : 371 bp, 106 bp
primer 1	primer 3	106 bp	KI	

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.