

RBRC 10186

Strain name C57BL/6-Phf6^{em1 Utr}

[illegible]

Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions


step	temp (°C)	time (sec)	note
------	--------------	---------------	------

[illegible]

Reagent	Conc	Reaction 1	2	3	4	
H ₂ O		2.8				(μL)
PCR buffer	2 X	10.0				
dNTP	2.0 mM	4.0				
primer 1	10 pmol/μl	1.0				
primer 2	10 pmol/μl	1.0				
DNA		1.0				
Taq	1 U/μl	0.2				
						Martini
						Gelatin
total		20				

Marker Gene Ladder 100 (Wako, Cat. 316-06951)

Gel	3% Nusieve in 0.5 X TBE (100 V, 50 min)
------------	---

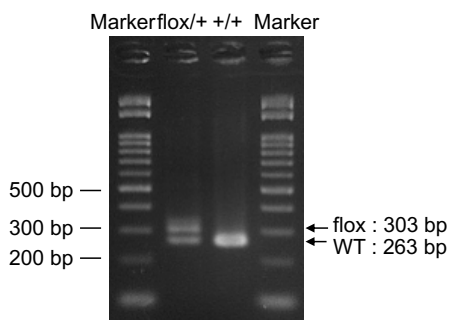

primer_set		product_size	product (allele type)	
primer 1	primer 2	263 bp	WT	+/+ (wild-type): 263 bp flox/+ (heterozygous): 263 bp and 303 bp
primer 1	primer 2	303 bp	flox	

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.